

Contact Information:

Ruchi Tiwari Genworks Health Pvt Ltd Unit 522-524, 5th Floor, Gamma Block Sigma Tech Park, Ramagondanahalli Varthur Road, Bangalore – 560066 Ph no. +91 96204 06263 ruchi.tiwari@genworkshealth.com

Ordering information

Item	Part Number
Echo-Screen Plus (TEOAE)	010127T
Echo-Screen Plus (DPOAE)	010127D
Echo-Screen Plus (AABR)	010127A
Echo-Screen Plus (TEOAE & AABR)	010127TA
Echo-Screen Plus (TEOAE & DPOAE)	010127TD
Echo-Screen Plus (DPOAE & AABR)	010127DA
Echo-Screen Plus (TEOAE, DPOAE, & AABR)	010127TDA
ES TS Set (TEOAE)	010127TS
ES DS Set (DPOAE)	010127DS
*Chaire of different learning and continued and affine death	harak Eska Carras arakan

Contact your local customer service representative for more information

Contact your local customer service representative for more in	ioiiidtioii.
Technical specifications	
Weight	Echo-Screen including battery pack and probe: 550g (1.2 lbs)
Size	Echo-Screen incl. battery pack: 230 x 95 x 53 mm (9.06 x 3.74 x 2.09 inches)
Battery	6V 1500mAh NiMH, exchangeable
Charger	100 - 240 V AC +/- 10% 50 - 60 Hz, 200 mA charge time approx. 2 hrs.
Power consumption	1.2 W max
Power management	auto backlight control, auto switch off, 2-level battery voltage sense, battery saver
Operating time	> 10 hrs with fully charged battery
Display	128 x 64 dot graphic LCD w/ switchable backlight H x W = 3.96 x 5.87 cm (1.56 x 2.31 inches)
Data memory	128 kB built-in flash memory, unlimited storage time
Interfaces	RS232 up to 115 kbps, infrared (optional modem feature available)
Storage temperature range:	0 – 40°C (32 – 104°F)
Operation temperature range:	5 – 35°C (40 – 95°F)
Humidity:	20-80 % rel.
Device class	lla (according to Council Directive 93/42/ EEC Appendix IX)
Applied standards	EN 60601-1 + A1 + A2 EN 60601-1-2 EN 60601-2-26 EN 60601-2-40
Note: Specifications are subject to change without notice.	



Robert-Koch-Str. 82152 Planegg

Tel: +49 (0)89 / 83942-0 Fax: +49 (0)89 / 83942-777

© 2012 Natus Medical Incorporated. All Rights Reserved. All product names appearing on this document are trademarks or registered trademarks owned, licensed to, promoted or distributed by Natus Medical Incorporated, its subsidiaries or affiliates.



Echo-Screen[®] **Hearing Screener**

The truly automated all-in-one hearing screener

Only the Echo-Screen Hearing Screener combines advanced TEOAE, DPOAE and AABR® technologies with fully automated results and completely standardized, pre-set screening parameters.

Fully automated results with fixed screening parameters

All results are based on the same standardized, objective criteria:

- Automatic pass/refer results that require no interpretation
- Preset screening parameters that cannot be adjusted Minimizes the risk of user error and test variability

TEOAE, DPOAE and AABR Technologies with patented binomial statistics algorithm

Transient Evoked Otoacoustic Emissions (TEOAE)

Trusted by clinicians for more than ten years:

- · Able to screen for even a mild hearing loss
- Testing frequency range: 1.5 to 3.5 kHz

Distortion Product Otoacoustic Emissions (DPOAE)

Frequency specific screening: 2, 2.5, 3.2, and 4 kHz

AABR Technology

- Necessary for detecting auditory neuropathy¹
- Optimized for screening newborns
- Testing frequency range: 2 to 4 kHz



KEY MESSAGES

- Hearing loss should be screened preferably before 1
- Universal neonatal screening rather than targeted 'high risk' screening is ideal.
- If abnormal OAE detected, it is repeated at 6 weeks or on the first immunization visit. If again abnormal, ABR
- is done for confirmation followed by full audiological evaluation and remediation with hearing aids.
- All NICU babies should undergo ABR testing to rule out auditory dys-synchrony/ auditory neuropathy



Robust signal detection

- Artifact rejection system distinguishes between external noise and actual patient response
- Fast test times
- TEOAE averages 20 seconds
- DPOAE averages 30 seconds
- AABR technology averages 2 minutes



AOAE Technology (TEOAE and DPOAE) tests a portion of the hearing pathway from the outer ear to the cochlea

COCHLEA

— = AABR[®] Technology

-- = AOAE® Technology

www.natus.com P/N 051745E www.natus.com

The **Echo-Screen** system utilizes proprietary AOAE® and AABR technologies

The Echo-Screen system meets all of your hearing screening requirements:

Perfect for any protocol and budget

- 17-key design provides convenient entry of patient information
- Multiple T (TEOAE), D (DPOAE), and A (AABR) technology combinations available
- Included Natus data management software allows data transfer
- 5-key design provides simple screening with no data entry
- TEOAE or DPOAE technology options are available
- Natus data management software allows transfer of results to a PC
- Ideal for one or two-step screening protocols
- Meets JCIH, AAP, and NCHAM requirements for surveillance screening of newborns and infants1



Echo-Screen Plus (17-key) configuration

Low refer rates

- Reduces the number of repeat tests
- Eliminates unnecessary follow-up and tracking

Fast, flexible and easy to use

- · Graphic display shows statistical result pattern
- Keypad allows fast data entry
- Ear probe is lightweight and easy to clean
- Battery is exchangeable provides approximately ten hours of testing
- Many language options are available

Fully automated results and screening parameters


- · Minimizes risk of user error and test variability
- · Reduces training time for new users

1 References on file at Natus





Automated TEOAE PASS Result

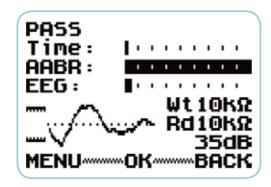

AABR Test Display



Printout of DPOAE Result



ECHO©SCREEN TDA + NoID 04/14 02:10 No Name 04/14 02:10pm

00,00,0000

DPOAE Result

06.05.04 12:34

(L2=50dBSPL)

2000Hz: PASS

2500Hz: PASS

3200Hz: PASS

4000Hz: PASS

Comment:

Examiner:

Right ear

OUT-Pat.





- Echo-Screen Plus (17-key) configuration:
- Stores up to 250 screening records
- Wireless printing of screening results
- Data transfer to and from a PC via Natus data management software
- Export compatibility with a variety of third party data management systems
- Echo-Screen Basic (5-key) configuration:
- Stores up to 50 screening records
- Printing of screening results via serial cable
- Transfer of screening results to a PC via Natus data management software

The Echo-Screen Plus system (17-key) includes*:

- Echo-Screen device
- Carrying case
- Ear tip accessory box
- Ear probe cable
- ABR cable (AABR configurations only)
- Two rechargeable batteries and charger
- Instructional materials (on CD)
- Natus data management software (on CD)
- *Y Probe cable and label printer purchased separately. Y Probe cable designed for use with Echo Couplers (AABR configurations only).

The Echo-Screen Basic system (5-key) includes*:

- Echo-Screen device
- Hard top carrying case
- Ear tip and accessory box
- Ear probe cable
- One rechargeable battery and charger
- Instructional materials (on CD)
- Natus data management software (on CD)

